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# Massage gun with articulated Joint User Manual

## INSTRUCTIONS FOR USE

Thank you for purchasing our SEECODE massage gun with articulated joint.

In order for you to enjoy your product for a long time and to use it to its full potential, please read these operating instructions carefully before using your device to avoid damage due to incorrect operation. Please pay particular attention to the safety information. If you give the unit to a third party, this instruction manual must also be handed over.

### Precautions

To prevent the risk of electric shock, fire and injury, this product must be used as follows.

1. do not reach for an appliance that has fallen into water.
2. do not place or drop into water or other liquids
3. never place or store the unit in wet environments or in high temperatures
4. do not drop the unit, use it improperly, remove the screws or attempt to disassemble, change or modify the unit
5. Do not operate the appliance for more than 30 minutes at a time.
6. always unplug the appliance from the mains before cleaning. Never operate or charge the unit unattended.

### Important safety instructions

1. the use of unauthorised voltage (DC 13V 1A) is prohibited.
2. use only suitable for adults, may only be used on dry body parts and clean skin. Please do not use over clothing.
3. use the massager only on soft tissues of the body, do not use on parts such as bones and head.
4. please only use the massage heads that are intended for the particular use.

5. strong pressure can cause bruising - please check the massage area constantly. You should stop using immediately if you experience any signs of pain or discomfort. Watch out for bruising.
6. Avoid contact with fingers and hair during use.
7. do not press hard on the massage head during use.
8. do not use the same attachment on the same spot for a long time.
9. infants, minors, pregnant women, persons with pacemakers, fractures or metallic parts in the body should not use the appliance.

10. persons with the following conditions should not use the massage gun:

- 1) Blood stasis
- 2) tumours
- 3) Heart disease e.g. with a pacemaker
- 4) If you have undergone joint surgery within 90 days, the product must not be used on this surgical site of the body and in the vicinity of about 10 cm


11. sensitive areas are: Head, face, cervical vertebrae, vertebral bodies of the spine, please be very careful here.

12. avoid positions close to the bone with little muscle coverage: e.g. the back of the foot and the back of the hand.

13. do not repair the massager yourself and do not replace any parts. Please contact customer service if you have any problems.

14. if a high temperature or foul odour should occur while using the appliance, stop using it immediately and unplug the power cord.

## Product description

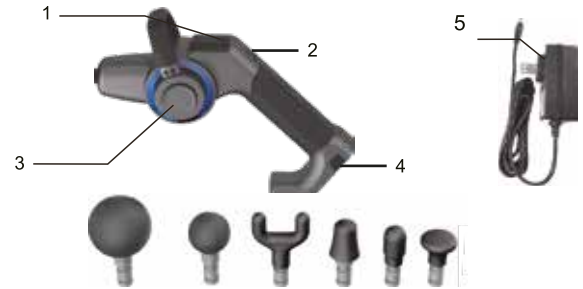
6 Massage-Heads

1. Speed control  
On / Off button
2. Charging port
3. Knob for angle adjustment
4. Compartment for batteries
5. Adapter

## Operating instructions

### Battery charge level and speed indicator

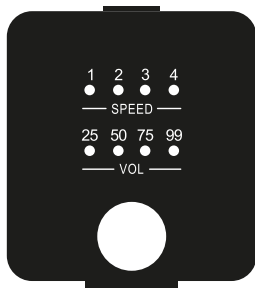VOL/ SPEED indicator

#### Before use

Fully charge the massage gun for approx. 2.5 hours before first use.

The massage gun runs for up to 4 hours on a full charge.

1. Plug the adapter into the charging port. The VOL indicator shows the battery level as follows:

4 red lights (=100%)

3 red lights, 1 red light flashing (<100%, ≥75%).

2 red lights, 1 red light flashing (<75%, ≥50%)

1 red light, 1 red light flashes (<50%, ≥25%)

1 red light flashing (<25%, ≥5%)

2. unplug the mains adapter from the massage gun and from the mains socket after it has been fully charged.

3. the battery must be fully recharged after three months of storage without use.

#### Notice:


The massage gun is only designed to operate with the rechargeable battery. It does not work when the mains adapter is connected.

You can also charge the batteries via the enclosed charging station so that you can continue to operate the massage gun with the second battery. To do this, connect the charging cable to the charging station and insert the battery into it from above. Make sure that the red dots on the battery and on the charging station point in the same direction. Only when the battery is firmly locked in after a noticeable resistance does the charging process begin.

### Using the massage gun

1. select a massage head and push it into the appliance with a twisting motion until it is fully inserted.

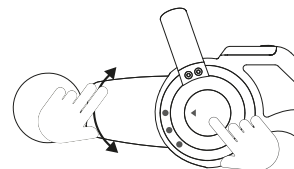To remove the head, pull it out of the appliance with a twisting motion.

2. press the angle adjustment button and select a desired angle of the massage head for yourself. There are 3 angles to choose from.

3. switch on the on / off button by pressing it a little longer and the battery charge level is indicated with red lights on the screen display. Note that the unit will switch off if it is not used for one minute, you will hear a warning tone. You must press the power button again to restart it.

4. Then press the on / off button again to select the speed in 4 levels. At the fifth press, the vibration turns off but the unit stays on.

5. when the massage gun is fully charged, the 4 LED lights are on, the 3 LED lights are on when the battery level is below 75%, 2 LED lights are on when the battery level is below 50%, 1 LED light is on when the battery level is below 25%.

When the battery level is below 10%, 1 indicator light flashes, please recharge the massage gun then.

6.To switch off the massager, press and hold the on/off button again. You can do this in any speed mode.

Note

A. The massager is equipped with an automatic switch-off after 15 minutes.

B. If you want to change the massage head, please always switch off the unit first.

C. Each time you press on the massage gun, an "indication tone" will sound.

D. During charging, the buttons do not respond.

### ***Massage head applications***

Small and large ball head

For large muscle groups such as arms and thighs

Air cushion head

For trigger points / pressure-sensitive hardenings in the tissue

Flat head

For all parts of the body

Fork head

For neck, spine and Achilles tendons

Ball head

For joint, deep tissue massage, precise point massage

### **Specifications**

Adapter:	Input: AC 100-240V 50/60Hz Output: DC 13V 1A
Power:	18W
Battery capacity:	2500mAh
Charging time:	approx. 2.5 hours
Operating time:	approx. 4 hours
Product size:	190 x 103 x 308 mm
Weight:	1kg

### **Package Contents:**

Massage gun  
2 rechargeable batteries  
6 massage heads  
AC adapter  
Storage box  
User manual

### **Troubleshooting guide**

#### **Problem**

The massage gun does not work

#### **Cause**

1. interrupted power connection or appliance plug causes poor contact. 2.

2. the power switch is switched off or the 15-minute timer is activated. 3. the motor has switched off automatically because if the appliance has been on too long: activation of the overheating timer.
3. the motor has switched off automatically because if the appliance has been on too long: activation of the overheating protection.

### **Solution**

1. Reconnect the appliance correctly to the power source.
2. switch the mains switch on again.
3. use the massager again only after the motor has cooled down naturally.

### **Problem**

The massage intensity is too high or too low.

### **Cause**

The massage head is positioned at the wrong angle.

### **Solution**

Adjust the angle of the head to reach an intensity that is moderate for you. Correct the speed mode to a weaker speed.

### **Guarantee/Warranty**

The warranty period is 24 months from the date of purchase. This includes a 6-month guarantee from the date of purchase. We grant a 6-month guarantee on the rechargeable battery.

The territorial scope of the warranty protection is Europe-wide. The proof of purchase is valid as proof of this. The product has been thoroughly checked for quality and functionality before dispatch. The guarantee is excluded in the following cases:

- in the event of modifications to the system without our approval
- in the case of use other than that stated above

- in the event of damage caused by falling

When making a claim under the guarantee/warranty, please return the complete product in transport-safe packaging to the point of sale from which the product was purchased. Describe the detected defect in the product and enclose the original proof of purchase with the shipment. The seller reserves the right to remedy the defect before refunding the purchase price. Necessary transport costs will be reimbursed to the customer if a justified claim under the statutory guarantee/warranty is established. The use of the legal rights for defects is free of charge.

CE Declaration of Conformity

Mobiset GmbH hereby declares that this product is in conformity with the essential requirements of the applicable European directives.

General information

Copyright

This document is protected by copyright. Any duplication or reprinting, including excerpts, as well as the reproduction of illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

Notes on environmental protection

The packaging materials used can be recycled. Dispose of packaging materials that are no longer needed in accordance with the locally applicable regulations. This product must not be disposed of in normal household waste within the European Union. Dispose of the appliance at municipal collection points.

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