**UV-C Dentalbox**  
**UV-C dental box**

**DE (2) | EN (8) | FR (14) | NL (20)**

## DE | Seecode UV-C Dentalbox

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## EN | Seecode UV-C dental box

Thank you for purchasing our SEECODE dental box, with which you can clean braces easily, quickly and hygienically. The UV-C disinfection effectively neutralises micro-organisms and bacteria. The rechargeable battery gives you complete flexibility in handling.

To ensure that you enjoy your product for a long time and can use it to the full, please read these operating instructions carefully before using the SEECODE dental box to avoid damage due to incorrect operation. Please pay special attention to the safety information. If you give the SEECODE dental box to third parties, these operating instructions should also be handed over. Please keep the operating instructions for future reference. Our product is not suitable for operation by children.

### **Scope of delivery**

Battery-operated dental box with removable inner part  
USB / USB-C cable  
Operating instructions

### **Technical data**

Dimensions: 13 x 9.2 x 4.2 cm

Weight: 108 grams

Battery type: Li-Ion 802540, 3.7V

Battery capacity: 800mAh


Charging time: up to 120 minutes, depending on battery charge level

Running time: up to 120 minutes


Material: ABS (scratch-resistant treatment) and silicone


## Product description

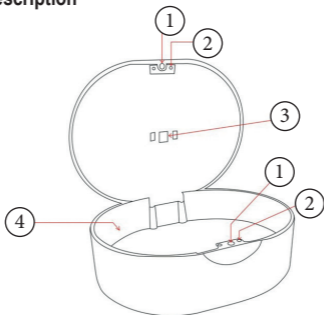1. Magnetic closure
2. Brass current contact sensor
3. UV-C lamps
4. Removable silicone inner part

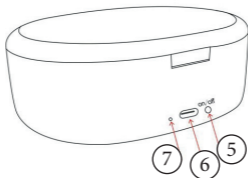5. ON/OFF switch
6. Type-C charging connection
7. Charge level and power indicator light

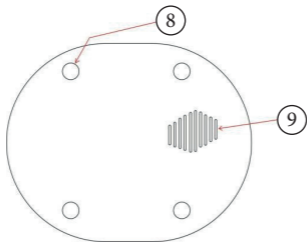- 8. Anti-slip pads
- 9. Fan opening

The material of the outer case is made of environmentally friendly and safe plastic, and the removable inner part is made of soft silicone to protect the braces from wear and tear.

The entire surface of the outer cover has a scratch-resistant seal to prevent damage during use.

The small brass contacts on the cover and box interrupt the cleaning process immediately when the cover is opened during use to prevent eye damage from ultraviolet radiation.

The ventilation is used for air exchange to keep the air inside the box dry. This is to reduce possible odors, e.g. from a worn retainer.

### **Operating instructions**

Place the brace in the silicone inner part of the dental box and switch it on with the small round button on the back. The UV-C light is then switched on and the ventilation starts to work.

When switching on and during cleaning, the operating indicator lights up green.

To switch off, press the button again and the green light goes out.

After switching on, the cleaning process is carried out for 15 minutes and then switched off automatically.

Note: For intensive cleaning, the dental box should be switched on several times in succession, each time at 15-minute intervals.

### **Charging the battery**

Plug the Type-C USB charging connector into the socket on the back of the unit and connect the other end of the cable to a 5V/2A power adapter (not included).

During charging, the charging indicator lights up red, when fully charged, the charging indicator lights up green.

Please fully charge the dental box before using it for the first time.

The charging time is up to 2 hours.

Please fully charge the dental box more often, at least once a month, to maintain continuous and good performance of the battery.

### **Safety warning for the built-in battery**

The battery may only be replaced by the manufacturer or an authorized service center.

Do not disassemble or open with tools.

Do not expose the dental box to direct heat or sunlight.

Keep the unit clean and dry.

If the charging port is dirty, please use only a dry cloth or brush to wipe it.

Please do not charge the unit steadily for a long period of time.

When using a temperature of  $20^{\circ}\text{C} \pm 5^{\circ}\text{C}$ , the battery has the optimal performance.

## **Care and maintenance**

The surface can be cleaned with a soft, if necessary slightly damp, cloth. To clean the silicone inner part, it should be removed. It can be rinsed under warm water and then reinserted after it is well dried. You can also rinse the brace well in the silicone inner part. Always dry the silicone inner part well before using it in the dental box.

Do not use abrasive cleaners or harsh chemicals for cleaning.

## **Warranty / Guarantee**

The warranty period is 24 months from the date of purchase. This includes a 6-month guarantee from the date of purchase. We grant a 6-month guarantee on the battery.

The geographical scope of the warranty protection is Europe-wide. The proof of purchase is valid as proof of this. The product has been thoroughly checked for quality and functionality before dispatch. The guarantee is excluded in the following cases:

- in the event of modifications to the system without our approval
- in the case of use other than that stated above
- in the event of damage caused by falling

When making a claim under the guarantee/warranty, please return the complete product in transport-safe packaging to the point of sale from which the product was purchased. Describe the detected defect in the product and enclose the original proof of purchase with the shipment. The seller reserves the right to remedy the defect before refunding the purchase price. Necessary transport costs will be reimbursed to the customer if a justified claim under the statutory guarantee/warranty is established. The use of the legal rights for defects is free of charge.

## CE Declaration of Conformity

Mobiset GmbH hereby declares that this product is in conformity with the essential requirements of the applicable European directives.

## General information

### Copyright

This document is protected by copyright. Any duplication or reprinting, including excerpts, as well as the reproduction of illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

## Notes on environmental protection

The packaging materials used can be recycled. Dispose of packaging materials that are no longer needed in accordance with the locally applicable regulations.

This product must not be disposed of in normal household waste within the European Union. Dispose of the appliance in accordance with regional regulations.

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Mobiset GmbH - Rösrather Straße 333 - 51107 Köln