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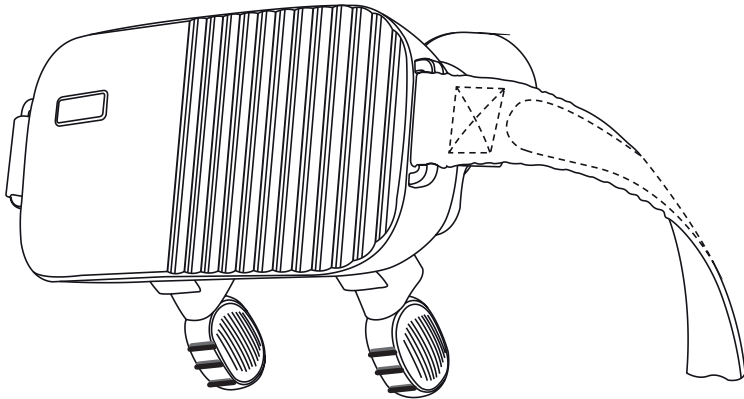
EN | SEECODE Neck and Shoulder Massager

Operating instructions

Thank you for purchasing our SEECODE neck and shoulder massager, which you can use to give your entire neck area a soothing and relaxing massage. As with a manual hand massage, 4 massage fingers press and stroke your neck and upper shoulder muscles.

To ensure that you can enjoy your SEECODE neck and shoulder massager for a long time and use it to its full potential, please read these operating instructions carefully and pay particular attention to the safety information before using the SEECODE neck and shoulder massager to avoid damage caused by incorrect operation.

If you give the SEECODE neck and shoulder massager to a third party, these operating instructions should also be handed over. Please also keep the operating instructions for future reference.

NOTE

The neck and shoulder massager is intended for use by adults only and must be kept away from minors.

Technical data

Battery: 2000 mAh lithium-ion battery

Charging cable: USB-C / USB-A connection cable, USB-C device connection (mains plug not included)

Power: 12 watts

Input voltage: 5 Volt / 2A

Heating levels 45, 50 and 55 °C

Remote control battery: CR2032 3.0V 210mAh

Size: 19 x 11.5 x 9.5 cm

Weight: 950 grams

Charging time: 4 hours when the battery is fully discharged

Running time: approx. 2 hours with uninterrupted use.

Scope of delivery

Neck and shoulder massager with straps

Battery charging cable

Remote control

Operating instructions


Important safety instructions, please note!

Keep the neck and shoulder massager out of the reach of children.

Do not rinse the device under running water and do not place it in water.

Do not use the neck and shoulder massager if it is defective in any way.

Productdescription

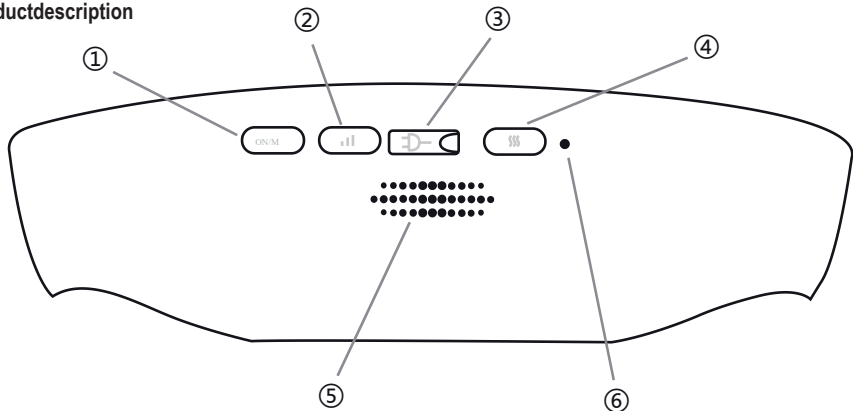1. On/off button and setting the modes
2. Button for setting the massage intensity
3. USB type-C charging connection
4. Button for switching on and setting the heat levels
5. Loudspeaker for voice announcements
6. LED charge level indicator
7. Signal lamp for heat function

Charging the battery

1. Connect the charging cable to the USB Type-C port on the device and a power adapter (not included)
2. Charging the neck and shoulder massager takes up to 4 hours, depending on the charge status of the battery.


The red LED on the charger indicates that the battery is charging. The green LED indicates that the battery is fully charged. The LED goes out when the device is disconnected from the mains.

Operating instructions

Place the massager on your neck and shoulder area.

Place the larger massage heads at the top of your neck and the smaller ones at the bottom, on the upper back area. You can press the massager onto your neck by pulling the straps forwards.

You can also use the strap to fasten the massager so that you do not have to use your hands to hold the massager. To do this, open the Velcro fasteners on the strap, guide the strap backwards under your arms and then cross it forwards to the abdominal area. You can secure the straps there using the Velcro fasteners.

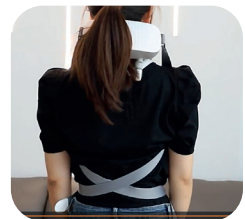Note: When switching the device on and off and for all settings on the device, the voice output confirms the respective mode.


A. Operating the massager itself

1. Press and hold the on/off button to start the massager.
The green light flashes when the device is switched on.
2. To start the device, please press the message intensity button (2).
Press the message intensity button to start the device at the lowest intensity level.
3. Briefly press the on/off button to switch to one of the preset modes:
Relaxing mode/ Awaken mode/ Combine mode/ Soothing mode.
Select the mode that is most comfortable for you and most relevant to the situation.
You can switch this mode off by briefly pressing it again.
4. You can also customise your own intensity level. To do this, briefly press the message intensity button and select the desired intensity by pressing it again:
Low intensity, medium intensity and high intensity.
Press briefly again to switch off the intensity mode.
5. For each mode, you can select 3 heat levels in the larger message heads by briefly pressing the heat button:
Low heat level / medium heat level / high heat level.
Press briefly again to switch off the heat mode.

B. Operation with the remote control

First remove the contact protection film from the battery by pulling it out of the remote control.

1. On/off button
2. Setting the message intensity
3. Setting the preset modes
4. Switching on and setting the heat function

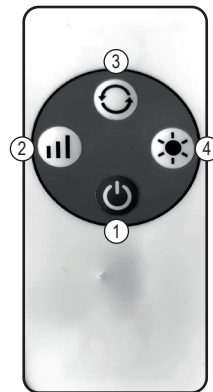1. **Before operating the remote control, press and hold the on/off button on the massager itself to start it.**
Then first press the message intensity button (2) on the remote control to start the device.
2. Use button 2 to select between low, medium and high intensity.
3. Button 3 is used to select between relaxation mode / wake-up mode / combination mode and calming mode.
4. Button 4 is used to select the low, medium and high heat level.
5. Switch off the device with button 1

Cleaning and care

Clean the device with a soft, slightly damp cloth if necessary. Do not use abrasive cleaning agents or harsh chemicals.

Do not use any metallic cleaning objects.

The massage heads can also be disinfected for hygienic cleaning.

General safety instructions for non-removable batteries

- Batteries may only be charged by children under adult supervision. Please keep away from children at all times.
- The electrolytes and electrolyte vapors contained in a rechargeable battery are harmful to health. Avoid direct contact at all times. If electrolytes come into contact with skin, eyes or other parts of the body, rinse immediately with plenty of fresh water. Then consult a doctor.
- Stop use or the charging process immediately if the battery swells, leaks, has a high temperature or emits smoke and strong odors.
- If you do not use the device for a longer period of time, charge it in between. Connected batteries can become deeply discharged and must then no longer be used or charged.
- Do not immerse batteries in water.
- Do not short-circuit the battery (connect + and - polarity).
- Do not disassemble or open batteries as this can lead to burns from the electrolytes.
- Do not expose batteries to excessive heat, e.g. direct sunlight, fire or similar. Do not place them on a heater.
- Do not dispose of rechargeable batteries by destroying or burning them and do not dispose of them in household waste. Dispose of them in accordance with local regulations.
- Please note that the end user is legally obliged to return used batteries and that used batteries can be returned to dealers free of charge.

Guarantee/warranty

The warranty period is 24 months from the date of purchase. This includes a 6-month guarantee from the date of purchase. We grant a 6-month guarantee on the battery.

The geographical scope of the warranty cover is Europe-wide. The proof of purchase serves as proof. The product has been thoroughly checked for quality and functionality before dispatch. The guarantee is excluded in the following cases

- if the product has been modified without our authorization
- use other than that specified above
- in the event of drop damage

If you make a claim under the guarantee/warranty, please return the complete product in transport-safe packaging to the point of sale from which the product was purchased. Describe the recognized defect in the product and enclose the original proof of purchase with the consignment. The seller reserves the right to subsequent fulfillment before reimbursement of the purchase price. Necessary transport costs will be reimbursed to the customer if a justified claim under the statutory guarantee/warranty is established. The utilization of statutory warranty rights is free of charge.

CE Declaration of Conformity

Hereby, Mobiset GmbH declares that this product is in conformity with the essential requirements of the applicable European directives.

General information

Copyright

This document is protected by copyright. Any duplication or reprinting, even in part, as well as the reproduction of the illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

Notes on environmental protection

The packaging materials used are recyclable. Dispose of packaging materials that are no longer required in accordance with local regulations.

This product must not be disposed of with normal household waste within the European Union. Dispose of the appliance in accordance with regional regulations.

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Mobiset GmbH - Rösrather Straße 333 - 51107 Köln